

HemoSeal

Tube Sealing System for perfect and reliable seals



High quality seals

Reliable operation

Quick and easy handling



Fresenius HemoCare

HemoSeal


The choice for controlled and safe tube sealing

HemoSeal is a desktop sealer, able to make continuous quality seals for PVC tubings up to 4,5 mm diameter. The sealing process is automatically triggered by placing a tube inside the head. Therefore operation is very efficient.

With an audio-visual alarm, the sealer is equipped with an reliable process control, that ensures leakage detection. The seals are well defined in shape and provide a rupture line, enabling the operator to separate the seals without contamination risk.



HemoSeal is provided with an additional adapter on the backside in order to use an extra hand sealer (hand sealer optional, order extra).

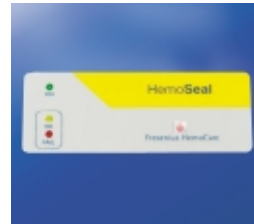

LED's for visual process control



The HemoSeal has a detachable head, which is very easy to clean.

Ordering information

For more information, literature, technical details etc. please contact your local sales representative.

HemoSeal

81039 HemoSeal International

Accessories

M638371 Hand Sealer

Technical data

HemoSeal	Dimension (WxDxH)	190x320x190 mm
	Weight	approx. 8 kg
	Power supply	110 V/230 V / 50 – 60 Hz
Hand sealer	Dimension	(WxDxH) 200x65x30 mm
	Weight	300 g
CE marked acc. to EMC- and LV-Directive protection against electrical shock:		class I

your local sales organization



Fresenius HemoCare

A company of the Fresenius Kabi AG

Fresenius HemoCare GmbH . Else-Kröner-Str. 1 . 61352 Bad Homburg . Germany
Tel. +49-61 72-6 08 8499 . Fax: +49-61 72-6 08 8498
www.fresenius-kabi.com